

Picovac 2<sup>™</sup> is a portable battery operated unit, able to create and maintain a selected pneumatic vacuum level for blood recovery systems.

The patient's own blood blood, normally lost during surgery, can be collected and reinfused to the same patient, thus avoiding the exposure to donated blood risks and diseases.



# **FEATURES** User friendly interface, with soft touch-waterproof keypad

## **Connections output:**

Female Luer Lock (male hose fitting 1/4" for D version)

### Possible pre-set instructions:

- Vacuum (mm Hg)
- Silence (no buzzer)
- Large backlighted graphic LCD with continuous readout of:
- Current vacuum value
- Set pressure value
- Use status
- Battery status

## Visual and acoustical warning for:

- Battery level
- Charging status
- Too low vacuum
- Too high vacuum



PICOVAC 2<sup>™</sup> bears the CE mark; in fact it complies with the requirements of Directive 93/42/EEC (as amended by 0476 2007/47/EC) for medical devices.


### **Lean Conformity Certificates.**


### **UNI EN ISO 9001:2015**

Design, manufacturing, distribution and service of electro-mechanical devices.

### **UNI EN ISO 13485:2012**

Design, manufacturing, distribution and service of electro-medical devices...

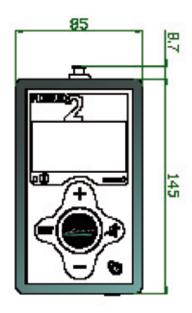

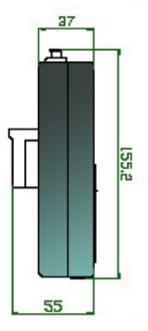

## **TECHNICAL DATA**

n° 2 AA pack type, Ni-Mh: 2,4 V 1300 mAh Batteries

at least 7 hours. Autonomy

not included, P/N: 0054700 (5 Vdc, 1.5A Medical Grade, **External Power Supply** 

Input 100-240V 50/60Hz)

Max Power consumption

Vacuum range  $-40 \div -110 \text{ mmHg}$  (0 ÷ -110 mmHg for D version)

Resolution 5 mmHq

Pressure measurement accuracy ±2 mmHg

(85 L x 160,5 W x 37 H) mm Dimensions

Net Weight 0,3 Kg

### Manufactured by:

### Lean s.r.l.

Via Marconi, 2 – 41036 Medolla (MO) – ITALY Tel. +39 (0)535 52200 - Fax +39 (0)535 51271

E-mail: lean@lean.it

